

BiomeOne®: multi-centric validation of a novel microbiome-based biomarker to predict response to cancer immunotherapy

Robinson I¹, Hochmair M¹, Schmidinger M², Ay L¹, Absenger G³, Pichler M³, Nguyen VA⁴, Richtig E⁵, Rainer B⁵, Jansen C⁶, Pacifico C⁶, Sladek B⁶, Knabl A⁶, Gasche N⁶, Valipour A¹

¹Department of Respiratory and Critical Care Medicine, Karl Landsteiner Institute for Lung Research and Pulmonary Oncology, Clinic Floridsdorf, Vienna, Austria; ²Department of Urology, Comprehensive Cancer Center, Medical University of Vienna, Vienna, Austria; ³Department of Internal Medicine, Division of Oncology, Medical University of Graz, Graz, Austria; ⁴Department of Dermatology, Medical University Innsbruck, Innsbruck, Austria; ⁵Department of Dermatology, Medical University of Graz, Graz, Austria; ⁶Biome Diagnostics GmbH, Vienna, Austria

INTRODUCTION

- The intestinal microbiome has a substantial influence on the cancer-related response to immune checkpoint inhibitor (ICI) therapy.
- BiomeOne® is a tumor-agnostic CE-IVD marked medical device that uses baseline stool samples to analyse the intestinal microbiome and predict patient response to ICI.
- The aim of our study was to evaluate the prognostic potential of BiomeOne® in a cohort of non-small cell lung cancer (NSCLC), melanoma and renal cell carcinoma (RCC) patients.

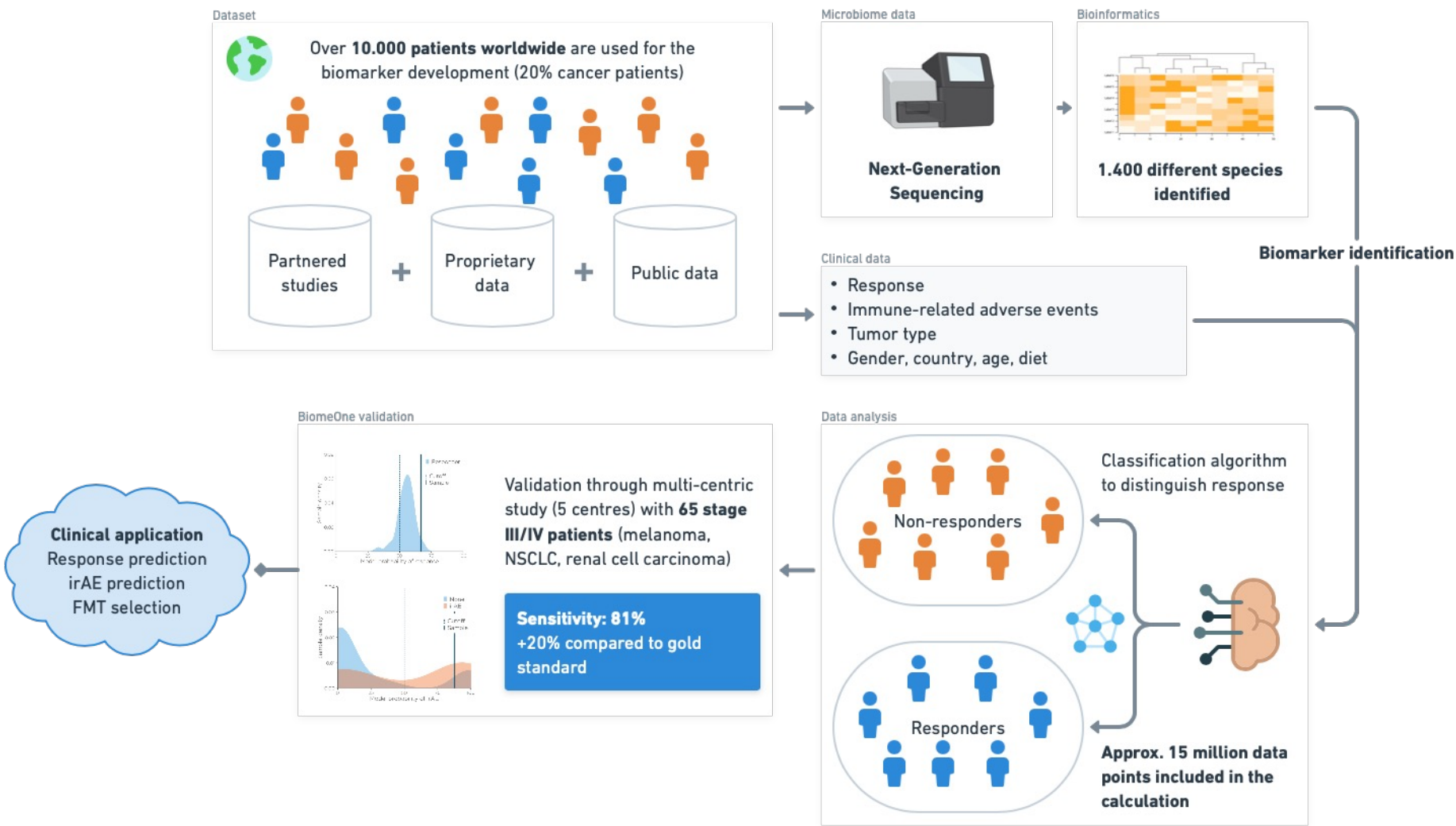


Figure 1. Schematic representation of biomarker development (=training dataset).

METHODS

- **65 patients:** NSCLC (42), RCC (16) and melanoma (7) treated with ICI. No systemic antibiotic treatment 30 days prior to ICI therapy initiation;
- **Stool samples** collected with at-home kit (Norgen Biotek), **prior to ICI treatment** initiation and at week 12;
- **16S rRNA sequencing** at both time-points;
- Statistical evaluation using Wilcoxon test, Mann-Whitney U test (alpha-diversity), PERMANOVA (beta-diversity) and ALDEx2, ANCOM-BC and MaAsLin2 (differential abundance);
- **Clinical response** was assessed at the end of first line therapy. **Responders (R)** were classified as complete and partial responders (CR, PR) and **non-responders (NR)** as stable and progressive disease (SD, PD). Best response was compared with the outcome of BiomeOne® analysis.

RESULTS

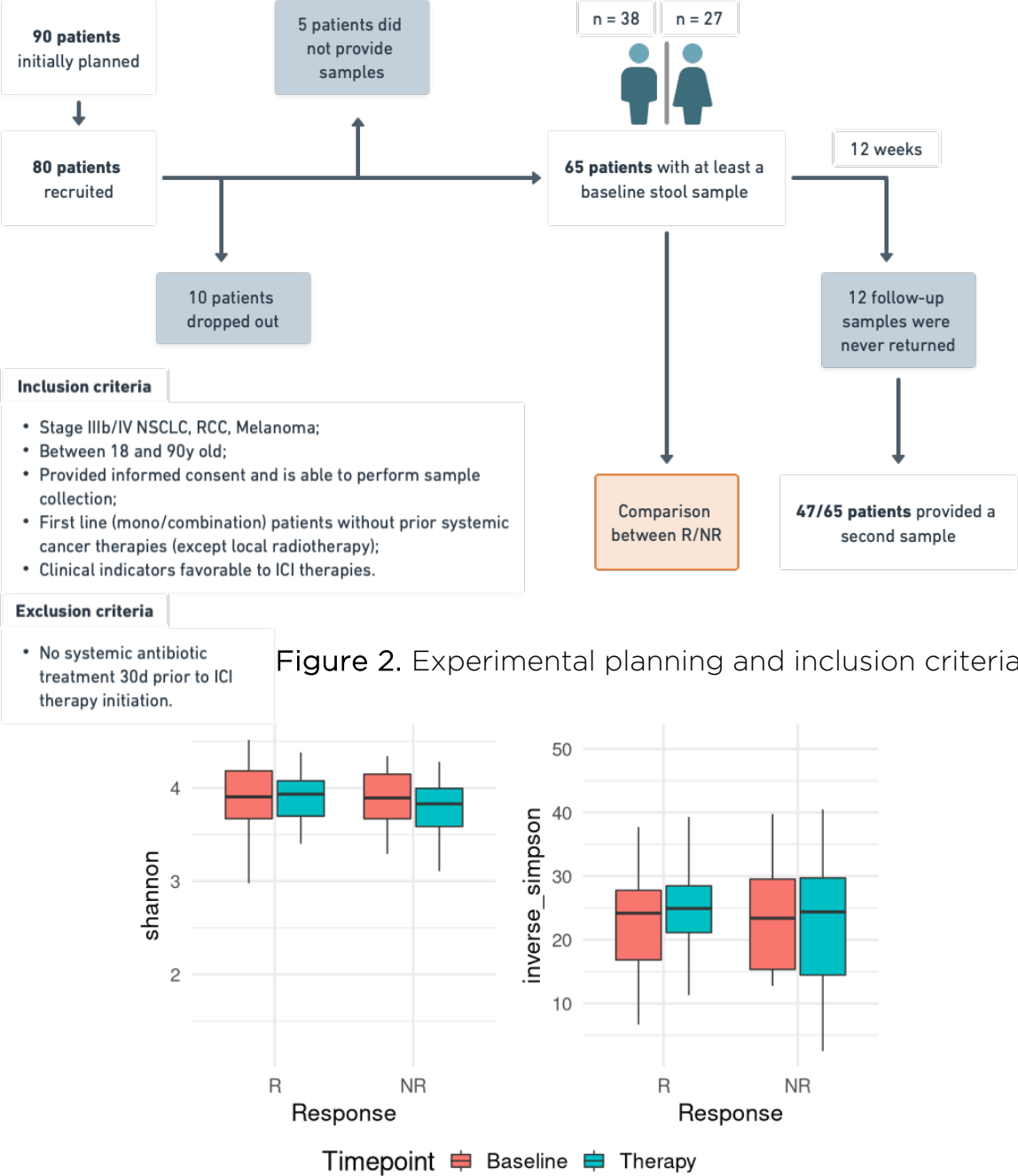


Figure 2. Experimental planning and inclusion criteria.

Figure 3. Alpha-diversity indices Shannon (P = 0.30) and inverse Simpson (P = 0.37) did not differ between R and NR nor between baseline and week 12 (P > 0.77).

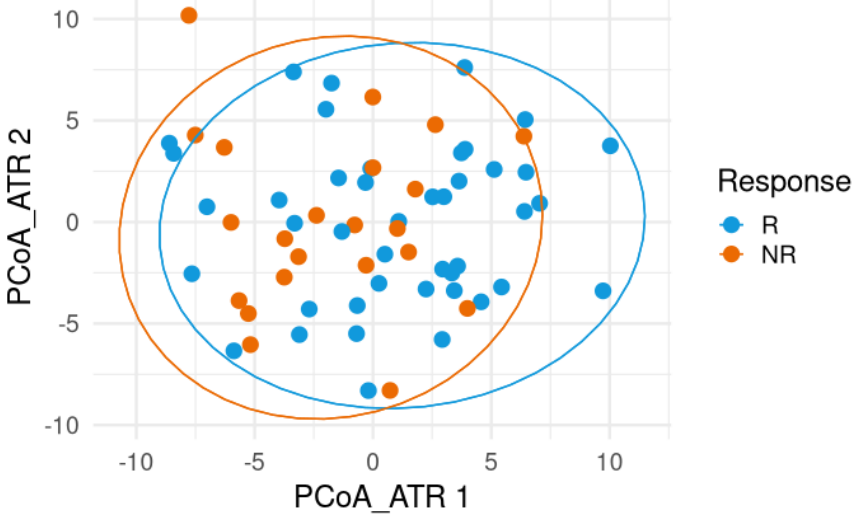


Figure 4. Beta-diversity did not differ between R and NR (P > 0.05).

- *Suterella* sp., *Faecalibacterium* sp. and *Haemophilus* sp. were found to be differentially abundant between R and NR (P < 0.05).
- No significant differences between the 3 tumour types, suggesting that the biomarker is tumour agnostic.

Table 1. Comparison between BiomeOne® results and clinical response.

	Responder CR + PR	Non-responder SD + PD	Total
BiomeOne: High	34 (pass)	10 (fail)	44
BiomeOne: Low	8 (fail)	11 (pass)	19
Total	42	21	63

- BiomeOne® was able to identify responders to ICI treatment with a **sensitivity of 81%** and a **positive predictive value of 77%**.

CONCLUSIONS

The presented research demonstrates the potential of BiomeOne® as a novel, non-invasive test to predict the outcome of ICI therapy. Further research will be necessary to expand the scope of BiomeOne® to other tumour entities.

Study sponsored by Biome Diagnostics GmbH. No conflicts of interest to declare.

Contact details
maximilian.hochmair@gesundheitsverbund.at

